



Panadent Articulating System



- Esthetics & Function
- Functional Denturism
- Diagnosis & Treatment Planning
- Smile Design & Facial Analysis
- Prosthodontics & Implantology
- Orthodontics & Oral Surgery

Product Catalog

Bite-Tab Compound Disks 180/box **Накусочные таблетки для фиксации прикусной вилки лицевой дуги.** **Features & Benefits**



- Constructed from red impression compound deposited on self adhesive strips.
- Can be used with most brands of bite forks.

Instructions:

- Remove Bite-Tabs by adhesive strip & place on bitefork.
- Temper compound in warm water before indexing on patients teeth.
- To remove used tab, dip in hot water and simply peel adhesive strip off of bite fork.

Bite-Tray регистраты 50/Bag

Features & Benefits

- Made of soft aluminum
- Designed for accurate jaw position interocclusal records
- Excellent for taking centric, protrusive and lateral border records.
- Can be used with any bite registration material



Универсальная линейка

Features & Benefits


The MMR is used for measuring tooth lengths, tooth widths, vertical overbite, horizontal overjet, as well as having opening and lateral range of motion scales.



Гайд для воскового моделирования по золотому сечению (set of 7)

Features & Benefits

Can be used with the Kois Adjustable Platform or by themselves to facilitate incorporating the golden proportion ratio of 1:1.618 for tooth widths into anterior restorations.

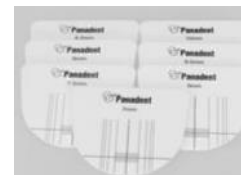

Беззубая челюсть без прикусного шаблона

Features & Benefits

The Panadent Edentulous Models are made of resin and include a bite record. They are used to simulate an edentulous patient and can be mounted to any articulator.



Поролоновый стабилизатор (60/пак)

Features & Benefits

The Bite Fork Stabilizers are used on the bottom of the bite fork so the lower teeth can support the bite fork against the upper teeth. It is a spongy material with an adhesive on one side for fixing to the bottom of the bite fork. Bite Fork Stabilizers can be used with any face bow system that utilizes a bite fork. It is clean and simple and eliminates the need for clumsy cotton rolls or messy bite registration material on bottom of bite fork.



Калибратор листовой пронумерованный

Features & Benefits

- Made of polyester strips
- Numbered (1-56) in 1/10th mm increments for convenience.
- Used as an anterior deprogrammer for taking centric records
- Can also be used in equilibration or coronaplasty procedures.



Линейка Золотого Сечения

Features & Benefits

The Golden Ruler measuring device is designed to help determine official dimensions, such as vertical dimensions, incisor edge position, and tooth width. It is calibrated to correspond to the "Golden principle (1:1.618) and is said to be useful for achieving aesthetic results.



Зубные ориентиры прозрачные (Set of 3)

Features & Benefits

The Clear Tooth Guides are made of clear acrylic blocks and can be used to evaluate tooth form, incisal edge and gingival contours.



Система денто-фациального анализатора Коиса

Features & Benefits

The Kois Dento-Facial Analyzer System registers and transfers the patient's occlusal plane as well as tilts in the occlusal plane in three planes of space to the articulator related to an average 100mm axis-incisal distance. O-rings can be adjusted on the Vertical Indicator Rod to the patient's facial landmarks for the evaluation of facial proportions.

The height of the lip commissures can be measured from the Index Tray for marking on the study casts.

Includes:

- #4305 Kois Dento-Facial Analyzer
- #4360 Kois Adjustable Platform



4300



4305

Пластины для анализатора Коиса (20/50/100)

Features & Benefits

The disposable index trays are used with the Kois Dento-Facial Analyzer and provides a permanent or disposable bite-fork registration record.



Столик для гипсовки анализатора Коиса

Features & Benefits

The Kois Adjustable Platform is used with the Kois Index Trays to mount the patient's maxillary cast to the articulator. The platform can be adjusted up or down to a new incisal edge or occlusal plane position.

Includes:

#4331 StandardWaxingGuide



Лицевая дуга Pana-Mount

The Pana-Mount™ Face-Bow features wrench or wrenchless detachable bite fork assemblies which allows the face-bow to be used for multiple patients without having to mount casts immediately. It has a standard 22mm nasion relator as well as an orbital pointer. Comfortable ear pieces and foam nasion pad increase patient comfort. The face-bow can be used as an average axis ear-bow as well as a true transverse axis bow (in conjunction with the Axis Mounting System and Axi-Path Recorder).

Includes:

Face-Bow, two complete Bite Fork Assemblies, Mounting Fixture, Standard Nasion Relator, T-Screws, and Hex Wrench



A

B

C

Bite Forks

Features & Benefits

Used with Pana-Mount FaceBow Constructed of all stainless steel, it is designed for easy cleaning and sterilization. The mid-line mark helps for orientation in the mouth.



Bite Fork Assembly (complete)

Features & Benefits

The Bite Fork Assembly is used with the Pana-Mount™ Face-Bow and comes complete with a standard bite fork and stem assembly, as well as T-Screws. By having extra bite fork assemblies, the operator can use one face-bow for multiple patients consecutively.



Magnetic PCH Articulator (#1610 AR)

Features & Benefits

The model PCH articulator is a precision articulator which allows the transfer of dental casts from one PCH or PSH articulator to another. The model PCH has a curved incisal pin calibrated in degrees. The curved pin holder allows the tip of the pin to stay in the same location on the incisal table when the vertical dimension is changed on the articulator. The positive centric latch depresses the centric pin to keep the articulator in perfect centric relation for mounting procedures or the pin can be depressed by hand to verify the centric position. The unique Dyna-Links keep the upper and lower frames joined together

movements.

Includes:

- 20 Mounting Plates
- Hex Wrench



PCF Magnetic Артикулятор

Features & Benefits

With an Integrated Magnetic System, Fixed Condylar Inclination set at 20 Degrees, Curved Incisal Pin, 1.5mm Analogs and Ability to Interchange Casts, the New Panadent PCF is a High Quality Articulator at a Great Price!

Includes: 10 Mounting Plates



PAL™ Plastic Articulator Lite (5/pcs)

Features & Benefits

- Used in conjunction with the Model PCH or PSH articulators.
- Creating a functional link between disposable and semiadjustable articulators.
- Has an anatomically correct axis of rotation; the same as the full size Panadent PCH or PSH articulator.
- Easy, lightweight and inexpensive alternative to multiple articulators.
- No super glue or molds required!
- Features include vertical and horizontal esthetic alignment guides as well as a built in, tilt back support stand.



Molded Plastic Case

Features & Benefits

This high impact, moulded Plastic Case is used to carry the articulator, face-bow, mandibular mounting stand, mounting plates and diagnostic casts with room left over for accessory supplies. Convenient to carry to the lab, study club or workshop meetings.

